



PELI
BIO THERMAL™

NanoCool™ SHUTTLE™ SYSTEMS



2-84141
Shuttle
Systems

COMPACT, LIGHTWEIGHT EVAPORATIVE COOLING SYSTEM

Shuttle™ Systems are refrigerated shippers that combine evaporative cooling technology for 2°C to 8°C overnight shipments of biological specimens and other healthcare products. Our evaporative cooling is seven times more powerful than gel packs or ice. It works by transferring internal heat to the outside through evaporation of a small amount of water stored at low pressure inside the package. Evaporation takes place at low temperatures and creates a powerful cooling effect. The active cooling process adjusts to variations in ambient temperature in a way that conventional cooling systems cannot.

Shuttle™ Systems have been qualified up to 42 hours and can be shipped inside FedEx or UPS Exempt Human Specimen (EHS) and UN3373 designated overnight envelopes. Shipping of UN3373 Biological Substance, Category B specimens inside the Shuttle requires the package and payload is tested in compliance with IATA Packing Instruction 650. EHS shipments do not require certification testing but do require marking the outer package or envelope with Exempt Human Specimen designation and following prescribed EHS packaging requirements.



DURATION
42 hours



TEMPERATURE
2°C to 8°C

WHY CHOOSE SHUTTLE™ SYSTEMS?

- Compact and lightweight
- Reduces transportation costs
- Qualified 2°C to 8°C for 42 hours
- Fits in overnight UN3373 designated envelopes
- Simple, one-touch activation
- No preconditioning hassles
- Uses evaporative technology

NanoCool™

PELIBIOTHERMAL.COM




SHUTTLE SYSTEMS PRODUCT SPECIFICATION

Product	Exterior Dimensions WxDxH "min" size due to angle		Usable Space WxDxH		Weight		Dim Weight
	mm	inches	mm	inches	kg	lbs	lbs
2-84121	305 x 200 (min 148) x 92	12 x 7.87 (min 5.82) x 3.62	190 x 90 x 30	7.48 x 3.54 x 1.18	1.09	2.40	3.3

SHUTTLE™ SYSTEM

SHUTTLE 42-HOUR PROTECTION OF 2-8°C PAYLOAD

Shuttle 42-Hours Chart
2-84121 Shipping System

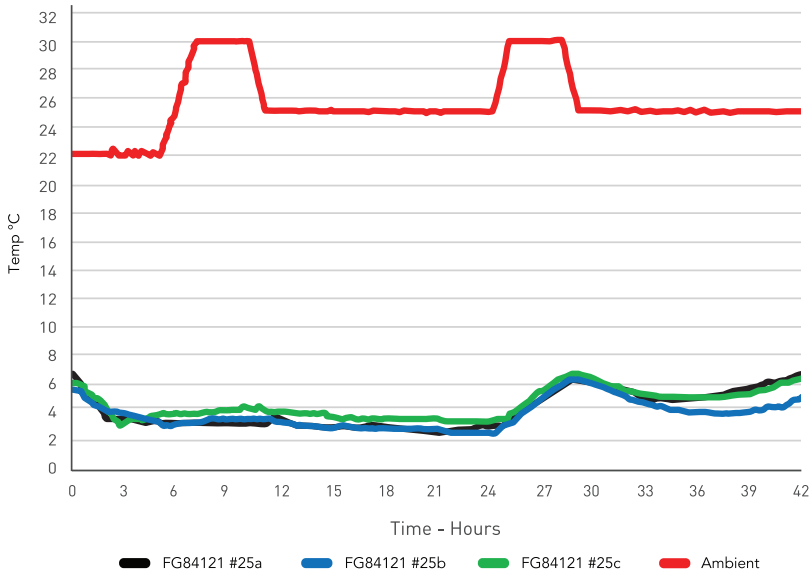

Qualification report available upon request

ONE-TOUCH ACTIVATION



NanoCool™

FOR MORE INFORMATION ABOUT PELI BIOTHERMAL PRODUCTS, PLEASE CONTACT:

EMEA OFFICE - INFONANOCOOL@PELIBIOTHERMAL.COM T: +44 (0)1525 243 770

AMERICAS OFFICE - INFONANOCOOL@PELIBIOTHERMAL.COM T: +1 484.716.9958

ASIA OFFICE - INFONANOCOOL@PELIBIOTHERMAL.COM T: +65 6681 0095

PELIBIOTHERMAL.COM

ISO17025, ISO9001-2008 Accredited

Golden Hour™, TIC™, Crêdo™, Crêdo Cube™, Crêdo ProEnvision™, Crêdo Xtreme™, Crêdo Duracube™, Crêdo ProMed™, Sherpa Systems™, CoolGuard™ and CoolPall™ are registered trademarks of Pelican Products Inc., its subsidiaries and / or affiliates



©2022 Peli BioThermal Ltd all rights reserved